

Project: ELISA paper project

Date: 20.05.2020

Performed by: Dorothy Smith

Approved by:

Test No.: 2

on: 22.05.2020 18:01

on:

Sample pipetting scheme

	Abbr	Description	Position
1	std1	Standard1 (0.0005 mg/mL)	C-4, C-5, C-6
2	std2	Standard2 (0.001 mg/mL)	F-4, F-5, F-6
3	std3	Standard3 (0.002 mg/mL)	A-7, A-8, A-9
4	std4	Standard4 (0.003 mg/mL)	C-7, C-8, C-9
5	std5	Standard5 (0.004 mg/mL)	F-7, F-8, F-9
6	std6	Standard6 (0.005 mg/mL)	D-10, D-11, D-12
7	std7	Standard7 (0.01 mg/mL)	
8	std8	Standard8 (0.015 mg/mL)	
9	std9	Standard9 (0.02 mg/mL)	
10	std10	Standard10 (0.025 mg/mL)	
11	std11	Standard11 (0.03 mg/mL)	
12	std12	Standard12 (0.04 mg/mL)	
13	std13	Standard13 (0.05 mg/mL)	
14	std14	Standard14 (0.06 mg/mL)	
15	std15	Standard15 (0.08 mg/mL)	
16	std16	Standard16 (0.1 mg/mL)	
17	std17	Standard17 (0.2 mg/mL)	
18	std18	Standard18 (0.3 mg/mL)	
19	std19	Standard19 (0.4 mg/mL)	
20	std20	Standard20 (0.5 mg/mL)	
21	std21	-	
22	std22	-	
23	std23	-	
24	std24	-	
25	std25	-	
26	std26	-	
27	std27	-	
28	std28	-	
29	std29	-	
30	std30	-	
31	sam1	-	A-1, A-2, A-3
32	sam2	-	B-7, B-8, B-9
33	sam25	-	
34	sam26	-	
35	sam27	-	
36	sam28	-	
37	sam29	-	
38	sam30	-	
39	sam31	-	
40	sam32	-	
41	sam33	-	
42	sam34	-	
43	sam35	-	
44	sam36	-	
45	sam37	-	
46	sam38	-	
47	sam39	-	
48	sam40	-	

Multiwell plate map

	1	2	3	4	5	6	7	8	9	10	11	12
A	BLQ	sam1	sam1	empty	empty	empty	std3	std3	std3	empty	empty	empty
B	sam4	sam4	sam4	empty	empty	empty	sam2	sam2	sam2	empty	empty	empty
C	empty	empty	empty	std1	std1	std1	std4	std4	std4	empty	empty	empty
D	empty	empty	empty	empty	empty	empty	sam3	sam3	sam3	std6	std6	std6
E	empty	empty	empty	empty	empty	empty	empty	empty	empty	empty	empty	empty
F	empty	empty	empty	std2	std2	std2	std5	std5	std5	empty	empty	empty
G	empty	empty	empty	empty	empty	empty	empty	empty	empty	empty	empty	empty
H	empty	empty	empty	empty	empty	empty	empty	empty	empty	empty	empty	BLQ

Initial measurement results

	1	2	3	4	5	6	7	8	9	10	11	12
A	BLQ	0.2682	0.2691	0.1939	0.1828	0.1705	0.6505	0.6937	0.6410	1.2601	1.2463	1.1666
B	0.0476	0.0495	0.0507	0.2441	0.2514	0.2441	0.7130	0.7221	0.6766	1.3088	1.2625	1.1298
C	0.0523	0.0509	0.0511	0.2877	0.2897	0.2910	0.8818	0.9079	0.9151	1.1787	1.2468	1.1355
D	0.0609	0.0574	0.0574	0.3544	0.3363	0.3568	1.0577	1.0606	1.0285	1.2611	1.3623	1.2885
E	0.0715	0.0691	0.0713	0.3796	0.3438	0.3722	1.0464	1.0070	1.0465	1.2791	1.3068	1.3249
F	0.0872	0.0833	0.0862	0.4594	0.4298	0.4366	1.1313	1.1287	1.0985	1.2950	1.3729	1.331
G	0.0997	0.1028	0.0918	0.4690	0.5093	0.4849	1.1799	1.1708	1.1049	1.2915	1.3207	1.3555
H	0.1168	0.1169	0.1163	0.5439	0.5892	0.5642	1.2311	1.1618	1.1438	1.2415	1.3131	BLQ

Calibration standards

Std. name	Number	Absorbance	Conc.	Variation
std1	3	0.2895	1.02	0.00136
std2	3	0.4419	1.04	0.01266
std3	3	0.6617	1.08	0.02293
std4	3	0.9016	1.22	0.01431
std5	3	1.1195	1.65	0.01489
std6	3	1.3040	7.39	0.04274

Model: Logit 5PL curve_fit fitting in lin_ln system

Absorbance = $D + ((A - D)/(1 + (\text{Conc}/C)^B))^E$

Model parameters

A= 0.149839, B= 1.610366, C= 0.030603, D= 1.564105, E= 0.252624

Model diagnostics

The Residual Sum of Squares **RSS** = 5e-05

Coefficient of Determination **R²** = 0.999616

Akaike Information Criterion **AIC** = -49.476479


Bayesian Information Criterion **BIC** = -50.517682

Coefficient of Correlation **r** = 0.999808


Time of calculatin = 0.043 [s]

Plots

Results for 5PL curve_fit in lin presentation

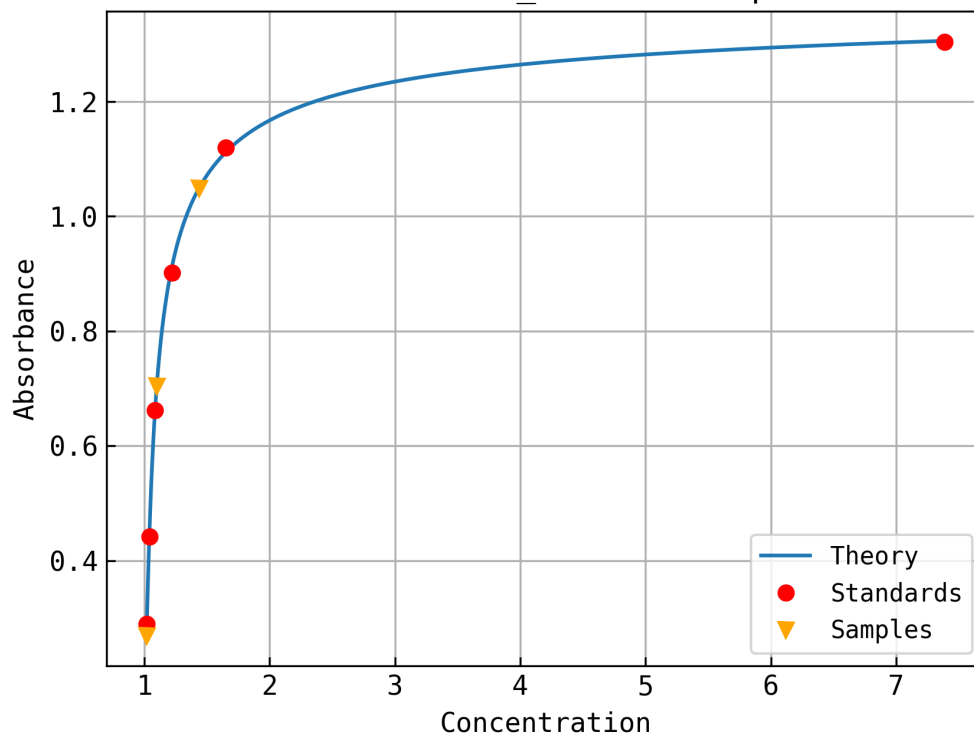


Results for 5PL curve_fit in lnx presentation

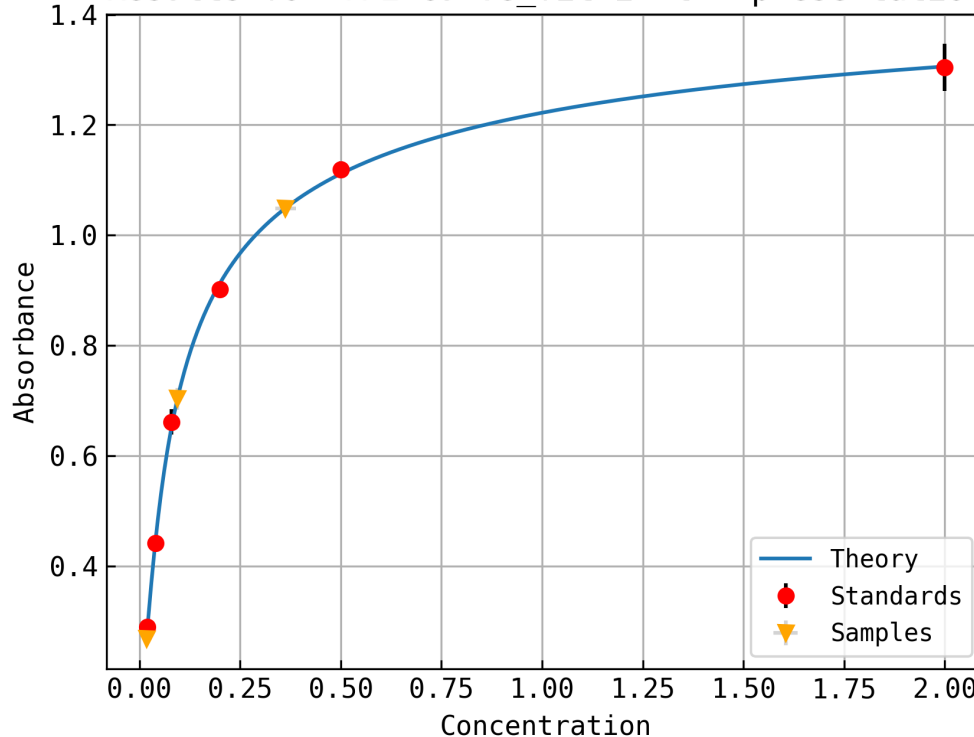


Plots

Results for 5PL curve_fit in lny presentation

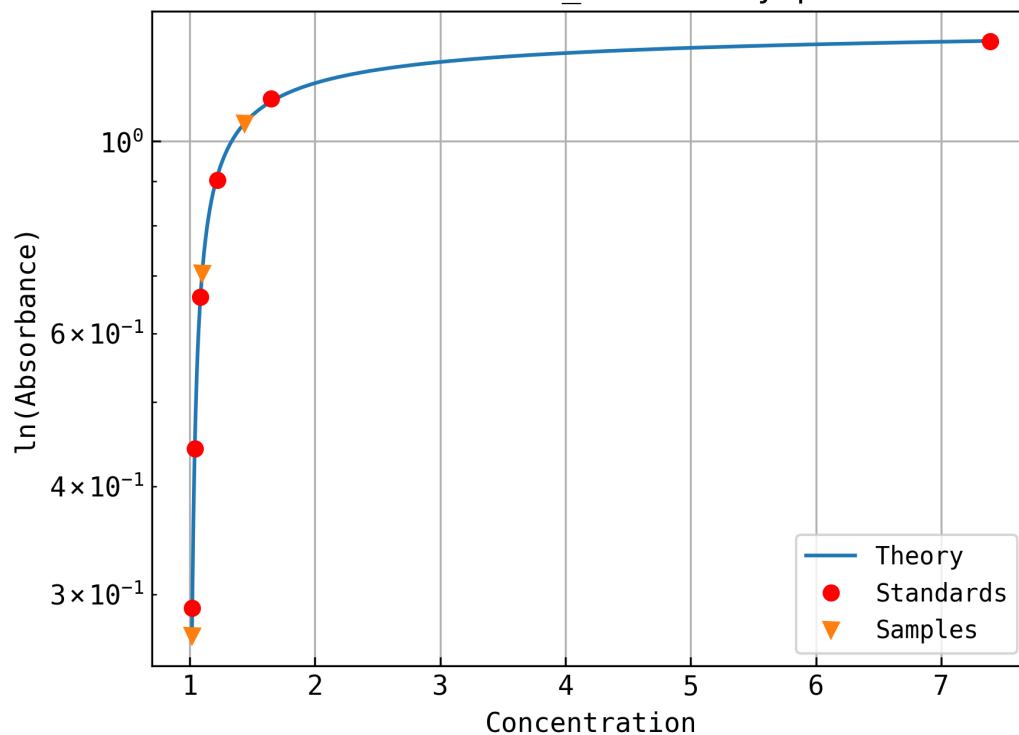

Results for 5PL curve_fit in lnln presentation

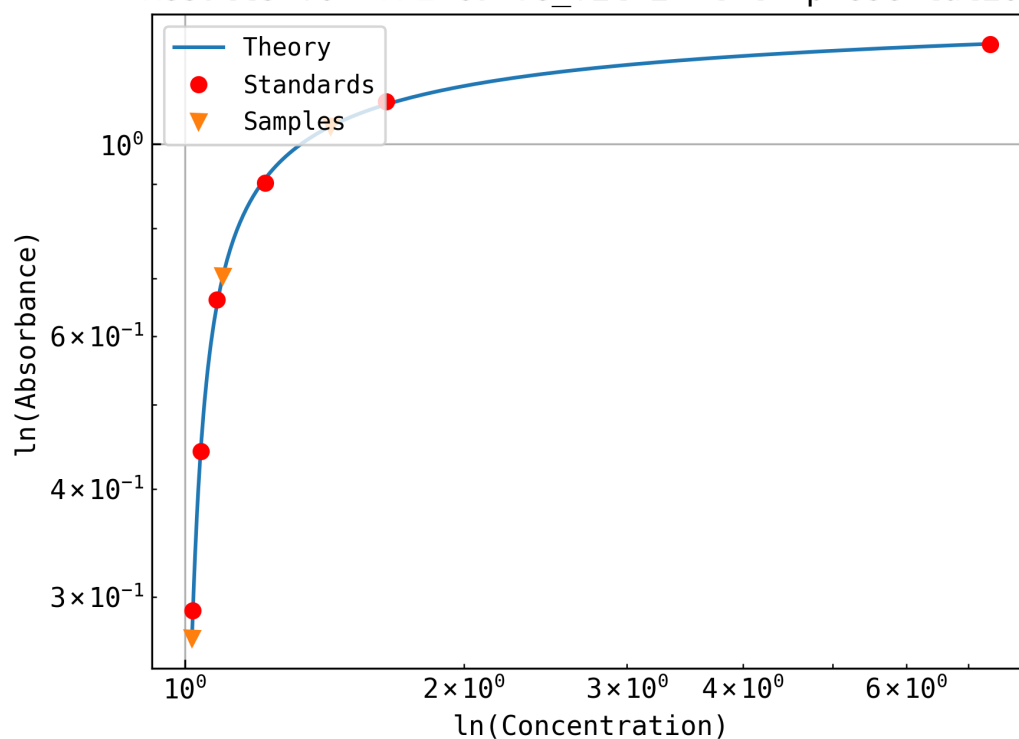

Calculation results

Good Samples	Conc.	SD down Conc.	SD up Conc.	Absorbance	SD Abs.
sam1	1.018	0.00005	0.00005	0.26865	0.00045
sam2	1.099	0.00641	0.00596	0.70390	0.01966
sam3	1.436	0.02680	0.02433	1.04893	0.01450

BLQ samples

	Bad Sample	Comment
1	sam1	blq_1_1
2	empty	blq_8_12
3	sam4	err1